



PetGuard Pro

Mobile Application User Manual

A Guide for Pet Owners

Prepared by the PetGuard Pro Team
(e21-3yp-PetGuard-Pro / MobileApp)

Contents

1. Introduction
2. Getting Started — Sign Up, Verify, and Log In
3. The Home Dashboard
4. Adding Your Pet's Profile
5. Analytics Tab
6. Features Tab
7. Location Monitoring
8. Health Monitoring
9. Activity Monitoring
10. Notifications
11. Profile & Settings
12. Help & Support
13. About
14. For Administrators
15. Troubleshooting Quick Reference

1. Introduction

PetGuard Pro is a real-time pet tracking and health monitoring mobile application. It pairs with a lightweight smart harness that streams live data to your phone, so you can:


- See your pet's live GPS location and set up safe zones (geofences)
- Monitor health vitals such as heart rate, body temperature, and respiratory rate
- Track daily activity, step counts, and detect falls or impacts
- Receive instant notifications when something needs your attention
- Generate PDF behaviour reports to share with your veterinarian

The app supports two kinds of accounts: Owners (pet parents who track and manage their own pets) and Administrators (who manage the overall user base and support tickets). This manual focuses on the Owner experience, with a short section on the Administrator dashboard near the end.

2. Getting Started — Sign Up, Verify, and Log In

2.1 Creating an Account

The first time you open PetGuard Pro you will land on the Login screen. If you don't have an account yet, tap Create new account at the bottom to open the Sign Up screen.

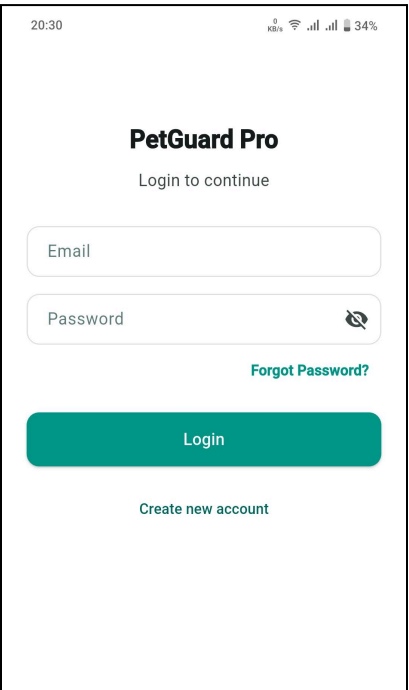

20:30 0 KB/s 34%

PetGuard Pro

Login to continue


Email

Password 

[Forgot Password?](#)

Login

Create new account

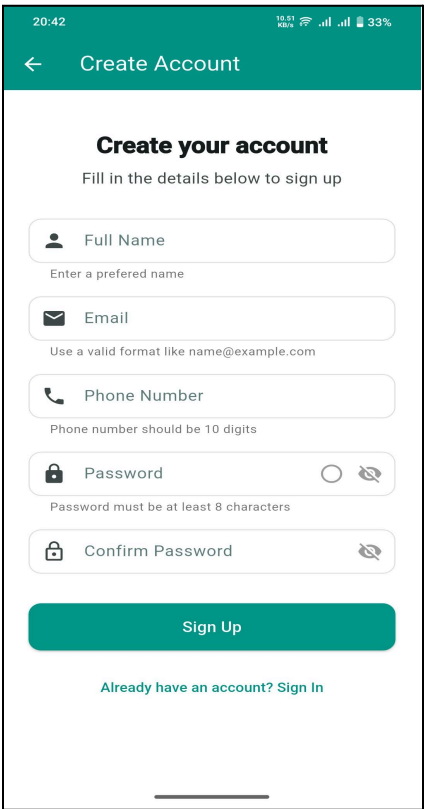

20:42 10.51 KB/s 33%

← Create Account

Create your account

Fill in the details below to sign up

 Full Name
Enter a preferred name

 Email
Use a valid format like name@example.com

 Phone Number
Phone number should be 10 digits

 Password 
Password must be at least 8 characters

 Confirm Password 

Sign Up

[Already have an account? Sign In](#)


fill in:

- Full Name — your preferred display name
- Email — a valid address such as name@example.com
- Phone Number — a 10-digit number (valid Phone Number)
- Password — at least 8 characters (a green check appears once this is satisfied)
- Confirm Password — must match the password above

Tap Sign Up. The app creates your account and immediately sends a verification link to your email.

2.2 Verifying Your Email

After signing up, you're taken to the Verify Email screen. This screen explains that a verification link has been emailed to you and will automatically continue once you click that link.

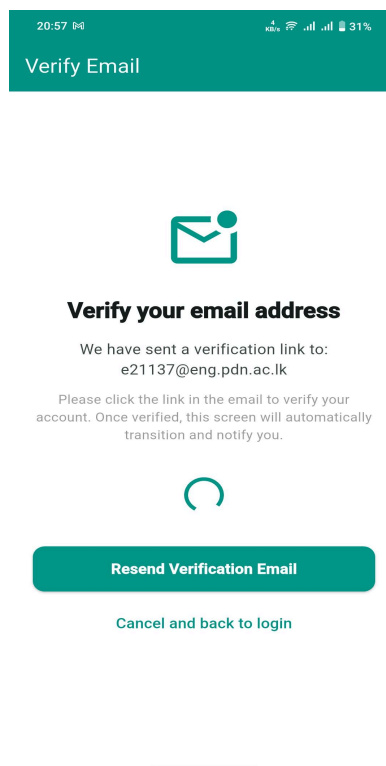

1. Open the email PetGuard Pro sent you.
2. Click the verification link inside it.
3. Return to the app — it detects verification automatically.
4. If you didn't receive the email, tap Resend Verification Email.


Note: after verifying your email, a PetGuard Pro administrator still needs to approve your account before you can log in and use harness features.

2.3 Logging In

Once your account is verified and approved, return to the Login screen, enter your Email and Password, and tap Login. The app remembers emails you've previously logged in with and offers them as autocomplete suggestions.

3. Adding Your Pet's Profile

In your first login to the application, the Home tab shows an "Add your pet details" instead of the vitals view. Tap Add Pet Details to open the pet profile form.

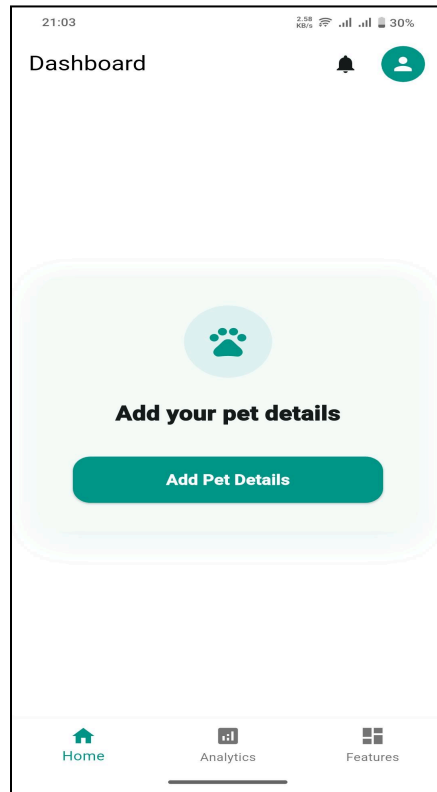


Fields on this form:

- Pet Name — a free-text name
- Size — Small (under 10 kg), Medium (10–25 kg), Large (25–45 kg), or Giant (over 45 kg)
- Age Group — Puppy (under 1 year), Adult (1–7 years), or Senior (over 7 years)
- Coat Type — Short/Hairless, Medium, or Long & Thick
- Flat Faced — Yes, No, or Not sure (relevant to breathing-related health thresholds)
- Activity Level — Low, Moderate, or High

Tap Add Pet Details to save. You'll see a confirmation message and be returned to the Home tab, which now displays your pet's live data. Now you have to wait until an admin assigns you a pet harness.

4. The Home Dashboard

The Home tab is the first screen you see after logging in. It's your at-a-glance summary of your pet's live status.

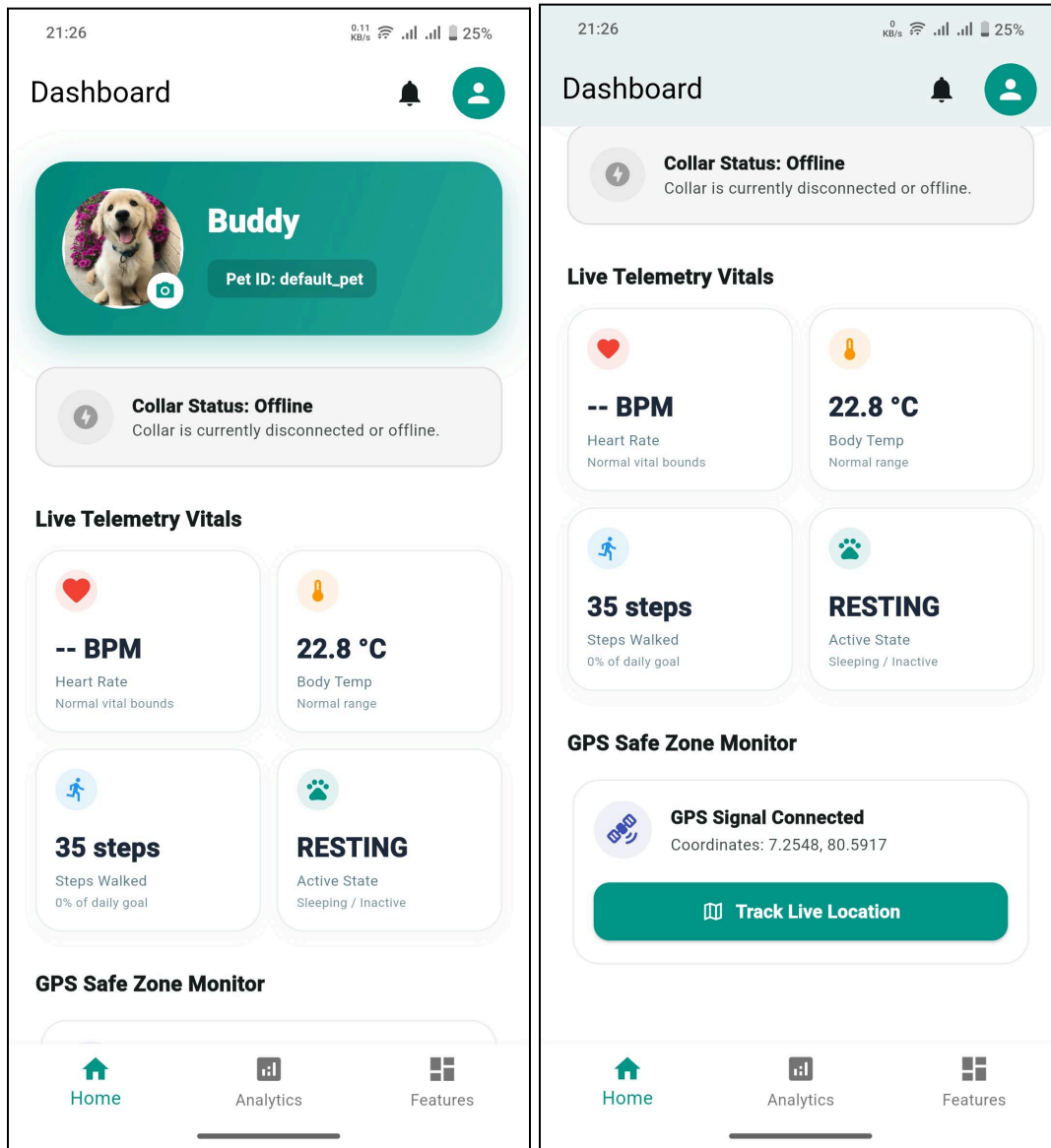

4.1 What's on this screen

- Pet card — your pet's photo, name, and breed. Tap the camera icon to upload or change the photo.
- Live Telemetry Vitals — four quick-glance cards: Heart Rate, Body Temp, Steps Walked, and Active State, each updated live from the harness.
- GPS Safe Zone Monitor — shows whether the harness has a GPS signal and its current coordinates, with a Track Live Location button that opens the location dashboard.

4.2 The top app bar

Every main screen shares the same app bar pattern:


-  Bell icon — opens Notifications, with a red badge showing the number of unread alerts.
-  Profile icon — opens a menu with My Profile, Settings, Help & Support, About, and Logout.

4.3 Bottom navigation

Three tabs are always available at the bottom of the Owner app: Home, Analytics, and Features. You can freely switch between them at any time — all live monitors (location, activity, health) keep running in the background regardless of which tab is open.

5. Analytics Tab

The Analytics tab (labelled Behavior Analytics) turns your pet's raw telemetry into a daily summary.

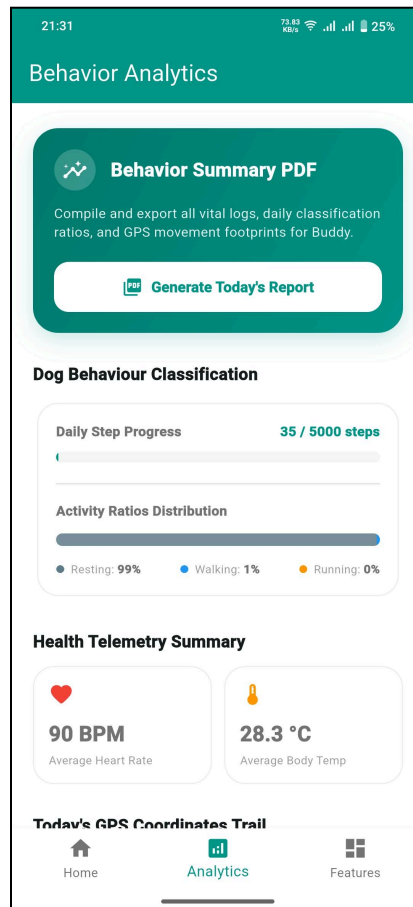

5.1 Behaviour Summary PDF

Tap Generate Today's Report to compile a PDF containing your pet's bio profile, daily behaviour classification, average health telemetry, and the day's GPS movement log. The file downloads to your device so you can save it or share it with your vet.

5.2 Dog Behaviour Classification


- Daily Step Progress — a progress bar comparing steps taken against a 5,000-step daily goal.
- Activity Ratios Distribution — a coloured bar showing the proportion of the day spent resting, walking, and running.

5.3 Average Vitals Today

Two summary cards show the day's average heart rate and average body temperature, each flagged if outside the normal range.

6. Features Tab

The Features tab (Real Time Features) is a simple launcher with three cards, one for each of PetGuard Pro's real-time monitoring systems. Tap any card to open its full dashboard.

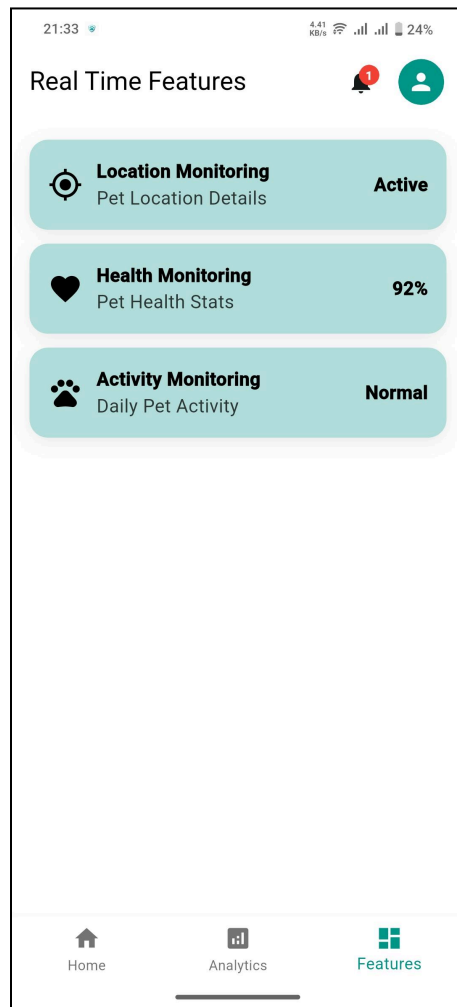


- Location Monitoring — opens the GPS Tracking dashboard (see Section 7)
- Health Monitoring — opens the Health dashboard (see Section 8)
- Activity Monitoring — opens the Activity Monitor (see Section 9)

7. Location Monitoring

The Location Monitoring feature covers everything related to where your pet is, has been, and the safe zones you've defined. It's the core feature of PetGuard Pro, giving you complete visibility and control over your pet's movements.

7.1 Location Monitoring Dashboard

This is your main hub for location monitoring. Access it from the Features tab's Location Monitoring card.

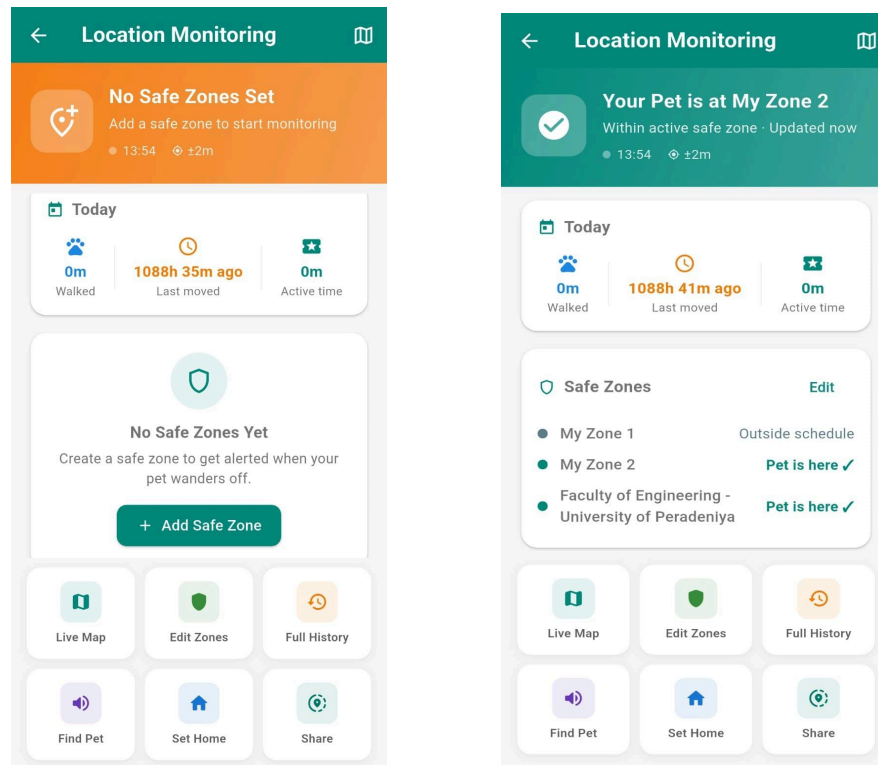

Status Banner (Top)

- Shows your pet's current safety status
- Colours indicate:
 - **Teal/Green** = Pet is safe (inside an active safe zone)
 - **Orange** = Pet is outside safe zones or in an off-schedule zone
 - **Red** = Pet has left a scheduled safe zone (alert triggered)
- Displays real-time accuracy ($\pm 2m$, $\pm 5m$, etc.) and timestamp of last update

Today's Statistics Card

Shows real-time tracking data for the current day:




- **Walked** (Paw icon): Total distance your pet has moved today
- **Last Moved** (Clock icon): Time elapsed since pet's last detected movement
- **Active Time** (Star icon): Total duration of active movement

Quick Action Buttons

1. **Live Map** (Map icon) – Opens full-screen interactive map with pet's current position.
2. **Edit Zones** (Shield icon) – Add, edit, or delete safe zones.
3. **Full History** (Clock icon) – Opens comprehensive location history.
4. **Find Pet** (Sound icon) – Remotely triggers the buzzer on your pet's collar.
5. **Set Home** (Home icon) – Creates an automatic ~100-metre circular safe zone, taking your pet's current GPS location as the center.
6. **Share** (Share icon) – Generates a shareable Google Maps link to your pet's current location.


Safe Zones Summary Card

- **Status:** Lists all configured safe zones – Shows zone name and status (e.g., "Pet is here ✓" or "Outside schedule")
- **Count:** Displays total number of safe zones created
- Tap the 'Edit' option to open the full Safe Zones management screen

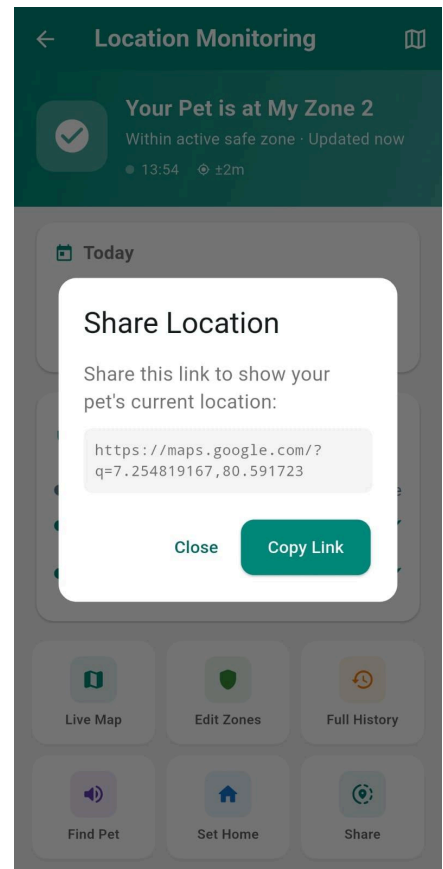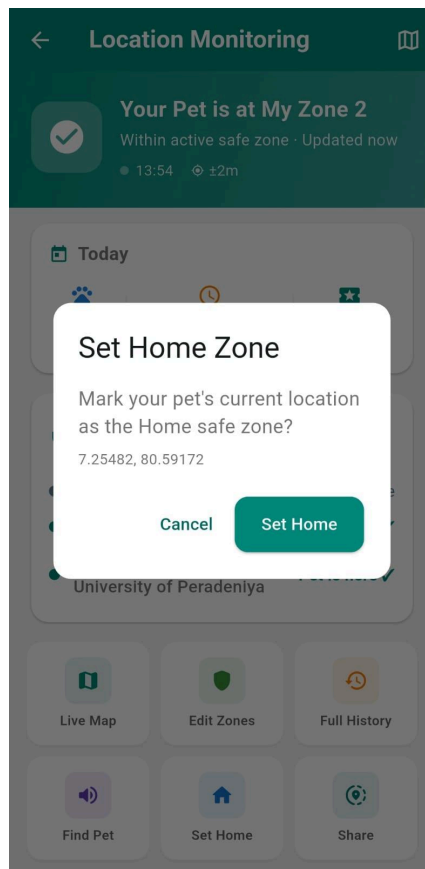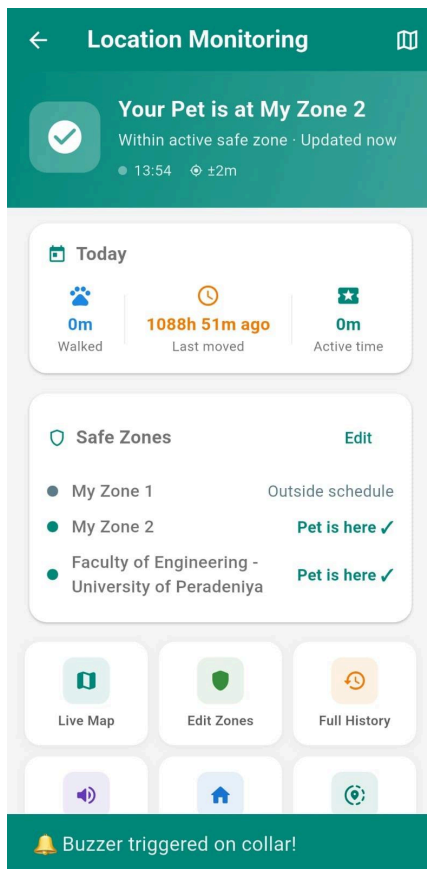

7.2 Live Map

Tap "Live Map" from the dashboard or map icon in the app bar to open the full-screen interactive map.

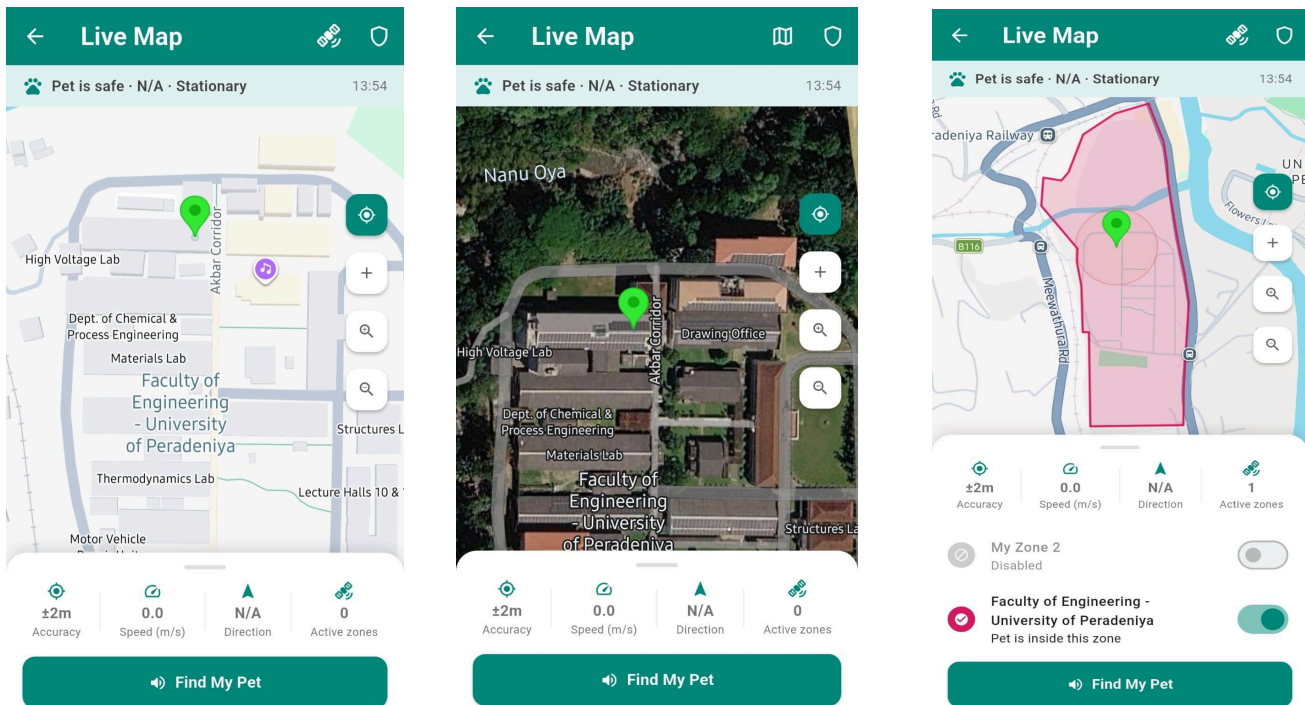

Map Display Features

Map Modes:

- **Standard Map** (default): Streets, labels, and landmarks visible
- **Satellite View**: Satellite imagery (toggle via satellite icon in top right)
- Works well for identifying specific buildings, parks, or geographic features

Pet Position:

- Large **green pin** shows your pet's current GPS location
- Centered on the pet's coordinates with the map zoomed appropriately
- Updates in real-time as GPS data arrives

Safe Zone Overlays:

- Each safe zone displays as a coloured boundary on the map
- Each zone labelled with its name
- Colour-coded for quick identification (Red, Green, Blue, Purple, Orange, Pink)

Pet Status Bar (Top):

- Shows current status: "Pet is safe · N/A · Stationary" (or moving)
- Displays time of last update
- Updates continuously as data comes in

Real-Time Data (Bottom Info Panel):

- **Accuracy ($\pm 2\text{m}$):** GPS accuracy radius in meters
- **Speed (0.0 m/s):** Current speed of movement
- **Direction (N/A):** Cardinal direction of movement (N, NE, E, SE, etc.)
- **Active Zones (0):** Number of zones currently active

Map Controls

Right-Side Control Buttons:


- **Center on Pet** (Crosshair/Target icon): Re-centers map on pet's current position
- **Zoom In (+):** Increases map zoom level
- **Zoom Out (-):** Decreases map zoom level

Find My Pet Button (Bottom):

- Large button spanning width of screen
- Tap to remotely trigger the collar's buzzer
- Useful when the pet is nearby but you can't see them

7.3 Managing Safe Zones (Geofences)

Safe zones (geofences) are virtual boundaries that alert you when your pet enters or leaves them. Create zones around your home, park, vet clinic, or any place you want to monitor.

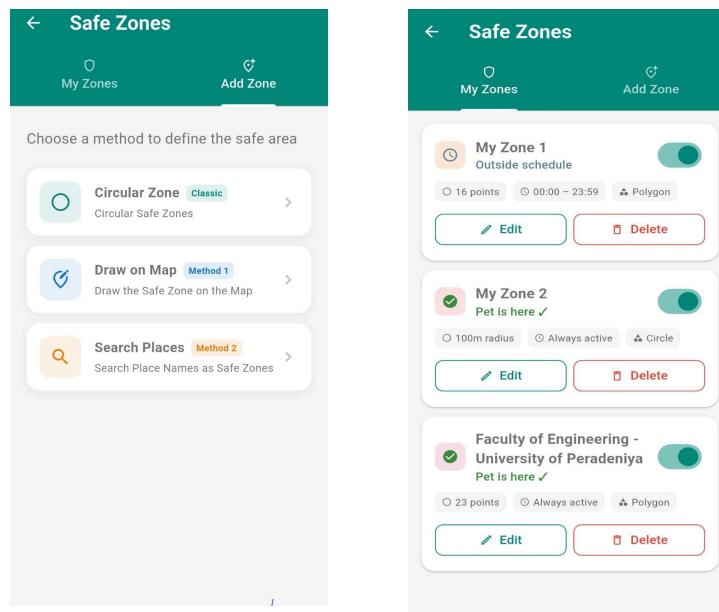

Managing Existing Zones

On Safe Zones List:

- Each zone shows a status card with:
 - Zone name and icon (colored circle)
 - Zone type: "Circle" (shows radius) or "Polygon" (shows number of points)
 - Status: "Pet is here ✓" or "Outside schedule"
 - Toggle switch to enable/disable
- **Edit** button: Modify zone settings, location, or size
- **Delete** button (red): Permanently remove the zone

Safe Zones Overview Screen

Tabs:



- **My Zones:** Lists all zones you've created
 - Shows zone name, size/type, status
 - Toggle zones on/off individually
 - Edit or delete existing zones
- **Add Zone:** Create a new safe zone

Creating a Safe Zone – Three Methods

Method 1: Circular Zone (Classic)

Best for: Home, parks, specific buildings


1. From "Add Zone" → Select "Circular Zone"
2. **Tap the map** to set the zone center point
3. A green pin appears at your tap location
4. **Drag the circle's edge** to adjust the radius
5. Typical radius: 10m-200m
6. Name your zone
7. Choose zone colour.
8. *Optional:* Set a schedule
9. Tap "Save Zone" to confirm

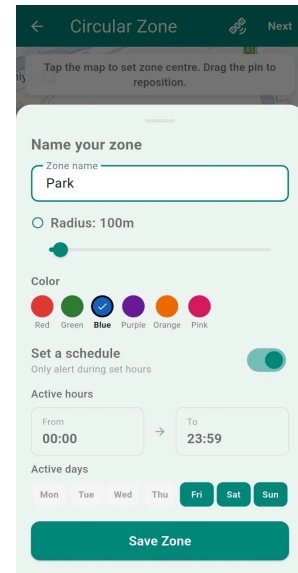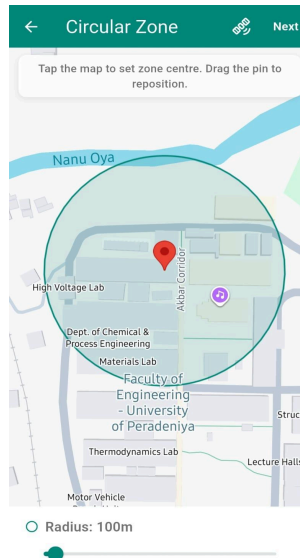

Method 2: Draw on Map (Polygon)

Best for: Irregularly shaped areas like yards, buildings, and neighbourhoods.


1. From "Add Zone" → Select "Draw on Map"
2. **Long-press the map** to place the first point (cyan pin appears)
3. Continue long-pressing to add more points around your desired area
4. Points connect with cyan lines showing your boundary
5. Connect back to the start point to close the polygon
6. When polygon is closed, it highlights in the zone color
7. **Adjust points** by dragging any cyan pin to reposition
8. Dialog appears: Name your zone, choose color, set optional schedule
9. Confirmation shows: "16 pts" (number of boundary points)
10. Tap "Done" to save the polygon zone

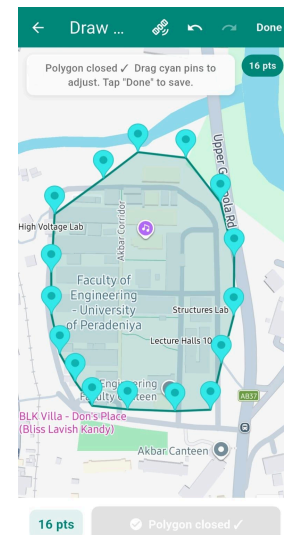

Method 3: Search Places

Best for: Named locations (restaurants, parks, landmarks, addresses)

1. From "Add Zone" → Select "Search Places"
2. **Type a location name** in the search bar
3. The app shows matching locations with red pins
4. Tap a location to select it
5. App displays the official boundary (teal outline) of the place
6. Shows: "Real boundary" badge confirming it's from Google Places data
7. Tap "Use this zone →" to create a safe zone using this exact boundary
8. Dialog appears to name the zone and choose settings
9. Tap "Save Zone" to confirm

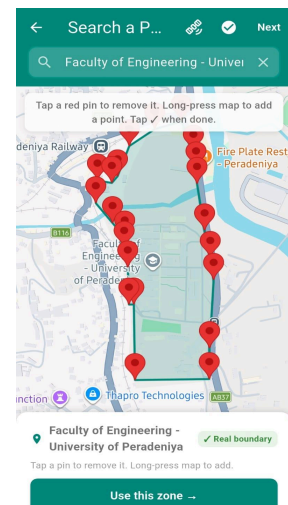


7.4 Location History

View a detailed history of where your pet has been. Access via the "Full History" button on the dashboard.

History Tabs


Tab 1: Map & Replay

- **Visual Route:** Shows a line connecting all GPS points for a selected time period
- **Red pins:** Mark individual location points
- **Pet trail:** Line shows the path your pet travelled
- **Replay Route** button: Animates the journey, showing movement over time
- **Time markers:** Show start time (13:45) and end time (13:54) of the route
- **Statistics shown:**
 - ❖ Total distance: "353m" (sum of all movements)
 - ❖ Average speed: "0.6 m/s"
 - ❖ Days covered: "1" (how many days of history shown)
 - ❖ Today's walked distance

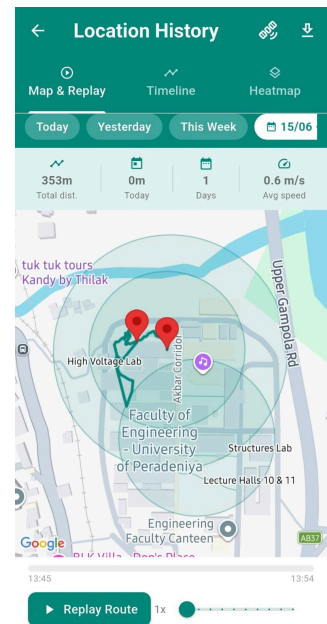


Tab 2: Timeline

- Chronological list of all location changes and zone events
- **Entries show:**
 - ❖ Time of event (13:45, 13:54, etc.)
 - ❖ Type of event (location update, zone entry, zone exit)
 - ❖ Zone name if relevant
 - ❖ For stationary periods, shows "last moved" timestamps
- Useful for understanding a pet's activity patterns throughout the day

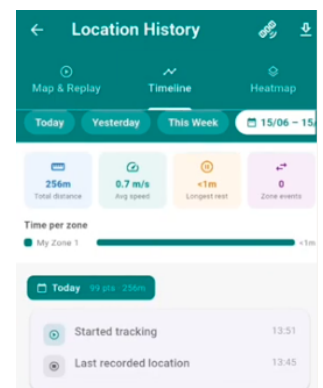

Tab 3: Heatmap

- **Colour-coded visualization** showing where the pet spends most time
- **Red areas:** High concentration (pet spends lots of time)
- **Pink areas:** Medium concentration
- **Light colours:** Low concentration (pet rarely visits)
- Overlays on the map to show activity patterns at a glance
- Useful for identifying favourite spots or unusual activity

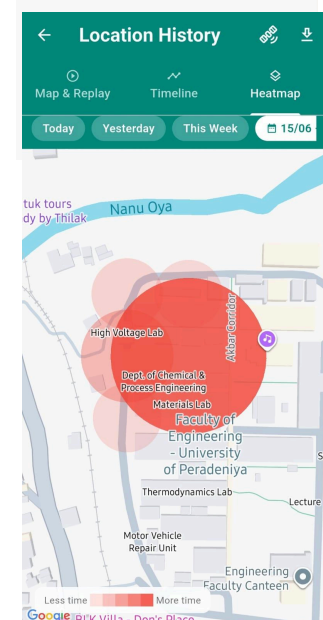

Time Period Selection

Filter buttons at top:

- **Today:** Shows location history for current day only
- **Yesterday:** Previous day's history
- **This Week:** Last 7 days of movement
- **Custom Date:** Tap calendar icon to select specific date range

Export Options


Export Location History Report dialogue (accessible from the history screen):

1. Spreadsheet CSV

- Exports data in comma-separated values format
- Opens in Excel or Google Sheets
- Includes: timestamps, coordinates, place names, zone information
- Contains: 538 location entries (example shown)
- Useful for analysis, record-keeping, or sharing with vet


2. PDF Report

- Professionally formatted report with:
 - Summary statistics (total distance, average speed, zones visited)
 - Place names and zone dwell times
 - Zone entry/exit timeline
 - Pie chart showing time distribution across zones
- Includes: location entries, place names, zone info, dwell time table
- Useful for vet visits, insurance claims, or professional records

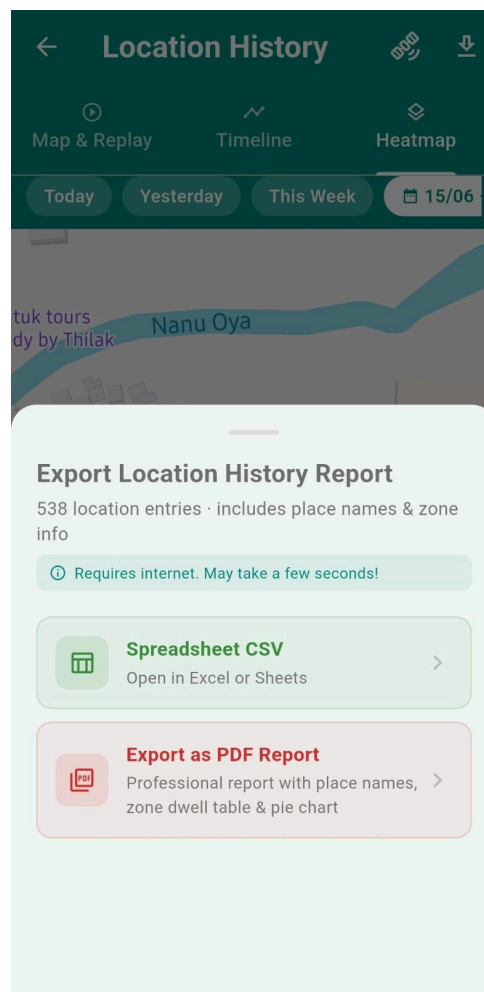

8. Health Monitoring

The health dashboard shows your pet's live vital signs, each compared against a normal range.

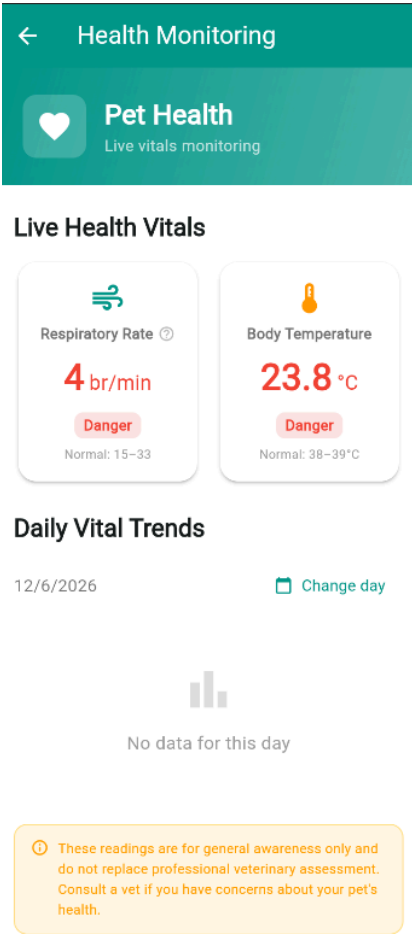

8.1 Vitals cards

- Respiratory Rate — breaths per minute; a key indicator of cardiovascular and lung health.
- Body Temperature — shown in °C, calibrated from the temperature sensor.

Each card is colour-coded — green for Normal, orange for Caution, and red for Danger — based on thresholds appropriate to your pet's profile. Tap the (i) info icon on a card for a plain-language explanation of what that vital means.


8.2 Trend chart

Below the vitals cards, a chart plots a selected vital sign over the course of the day. Use Change day to review a previous day's readings.

Alerts: If a vital crosses into the Danger range, PetGuard Pro sends you an immediate notification — see Section 10.

9. Activity Monitoring

The Activity Monitor tracks movement, steps, and sudden impacts (such as falls) using the harness's built-in motion sensor.

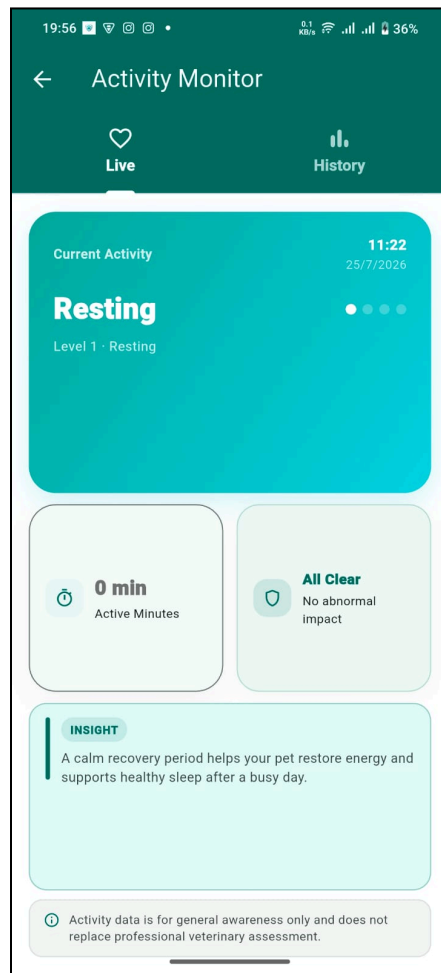

9.1 Live tab

- An impact status banner — green when nothing unusual has been detected, or a warning when an impact has been registered.
- Steps Today and Impacts Today summary cards.
- A 24-hour activity timeline chart, broken into blocks from 12 am to 11 pm.

9.2 History tab


Switch to the History tab to review past days' activity levels and impact events, useful for spotting behaviour changes over time.

What is an impact alert? The harness's G-sensor detects sudden high-acceleration events such as falls, jumps, or collisions. If an alert shows as triggered, check on your pet right away.

Live

History

 Today

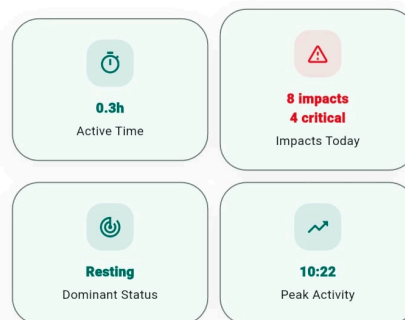


10. Notifications

Tap the 🔔 bell icon in the app bar from any main screen to see all recent alerts in one place — geofence zone crossings, activity impact detections, and health vital warnings. The red badge on the bell always reflects your current unread count and updates instantly as new alerts arrive — even if you're on a different tab.

11. Profile & Settings

Reach both pages from the 👤 profile menu in the app bar.


11.1 My Profile

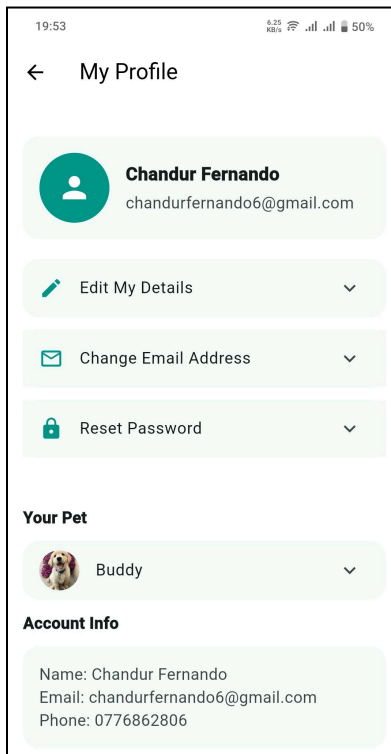

This screen shows your name and email at the top, followed by four expandable sections:

- Edit My Details — update your name and phone number
- Change Email Address — update your login email
- Reset Password — send yourself a password-reset email
- Edit Pet Details — update your pet's name, size, age group, coat type, and activity level

11.2 Settings

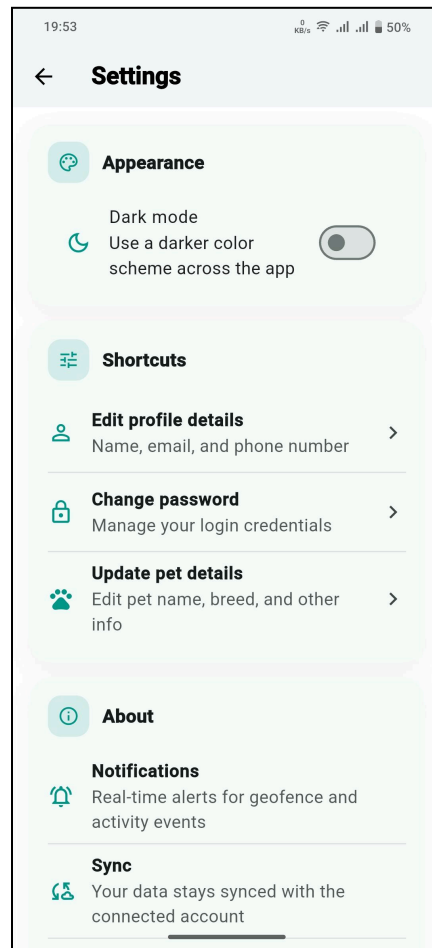

- Appearance — toggle Dark mode for a darker colour scheme across the app.
- Shortcuts — quick links to edit profile details, change your password, or update pet details.
- About — read-only notes on Notifications, Sync, and Privacy behaviour.

12. Help & Support

12.1 Frequently Asked Questions

Question	Answer
How do I configure geofence safe zones?	Go to Features → Location Monitoring → Edit Zones → Add Zone. Choose a circular zone or draw a custom polygon boundary.
How often are telemetry values updated?	The harness updates GPS, heart rate, and activity state every few seconds while your pet is moving. When stationary, it switches to a battery-saving mode with less frequent updates.
What are normal health vital thresholds?	Typical heart rate is 80–140 BPM for adult dogs and cats; body temperature typically ranges 37.5–39.2°C. You'll get an instant alert if a vital crosses outside these bounds.
What does the G-sensor alert mean?	The harness's motion sensor detects high-acceleration impacts such as falls, jumps, or collisions. If it shows as triggered, check on your pet immediately.

12.2 Troubleshooting Guide

Harness shows Offline status

5. Ensure the harness battery is charged (the LED should turn green).
6. Verify your home Wi-Fi network or cellular gateway is active.

No GPS location updates

7. Make sure the harness has an unobstructed view of the sky.
8. Allow up to 2 minutes for GPS to initialise outdoors from a cold start.

12.3 Submitting a Support Ticket

If the FAQs and troubleshooting steps don't resolve your issue, scroll to the bottom of the Help & Support screen and fill in the Submit Support Ticket form with your name and a description of the problem. Our support team will follow up by email.

13. About

The About screen (reachable from the profile menu) gives a short description of PetGuard Pro: an advanced real-time pet tracking and health monitoring solution that combines a lightweight smart harness with real-time cloud data to monitor location, log activity, configure geofence safe zones, and analyse vital-sign thresholds. It also lists the app's current version number and key feature highlights.

14. For Administrators

Accounts with the Administrator role land on a separate Admin Dashboard after login, with its own three-tab bottom navigation: Home, Users, and Support.

- Home — an overview of total users, registered pets, and open support tickets.
- Users — manage the full list of registered owners, view each pet's details, and assign a harness (pet ID) to a user's pet.
- Support — review, respond to, and close support tickets submitted from the Help & Support screen.

15. Troubleshooting Quick Reference

Issue	What to try
Can't log in	Confirm your email is verified and your account has been approved by an administrator. Use Reset Password from the Login screen's account flow if needed.
Harness shows offline	Check the battery/LED, confirm your network connection, then restart the harness (hold power for 3 seconds).
No GPS updates	Make sure the harness can see the sky; allow up to 2 minutes for a cold-start GPS fix outdoors.
Missed a geofence alert	Confirm notifications are enabled for PetGuard Pro in your phone's system settings and that the zone is still active under Edit Zones.
Vitals show "--".	The harness hasn't reported a fresh reading yet — this clears automatically once new telemetry arrives.

For anything not covered here, use the Submit Support Ticket form in Help & Support (Section 12.3), and our team will assist you.